

44

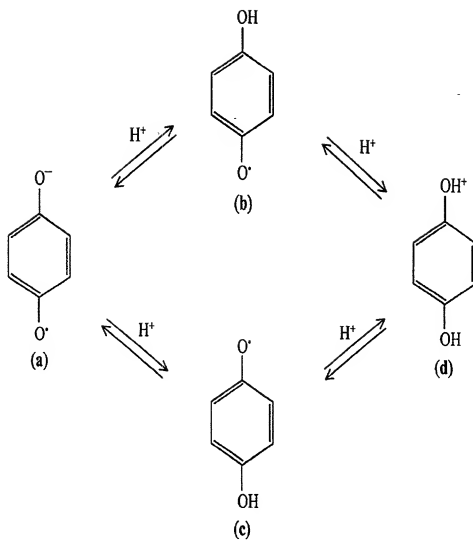


Fig 1. Mesomeric forms of p-benzosemiquinone. (a) anionic; (b),(c) neutral; (d) cationic.

10076033.052802

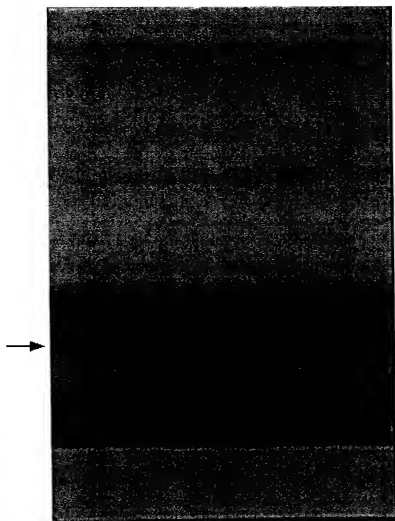


Fig 2. Band thin layer chromatography of the methanol solution after lyophilization (step 5). —> Indicates the band of cs-oxidant.

208250.5597001

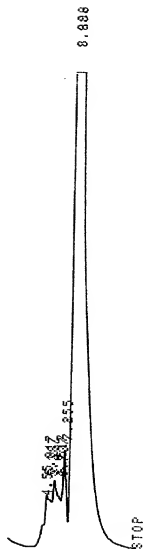


Fig 3. HPLC profile the butanol extract after TLC. The cs-oxidant (step 6) eluted as a major peak at the retention time of 8.808 min. The amount of cs-oxidant eluted was  $\approx 12 \mu\text{g}$ .

208250.6697001

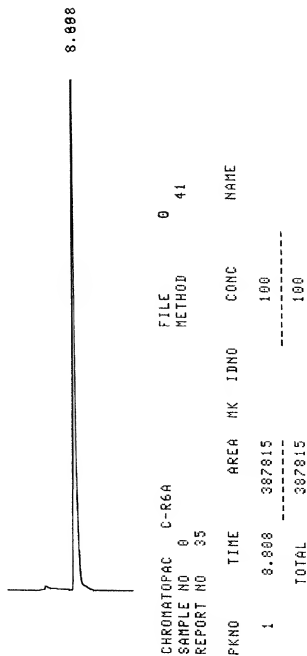


Fig 4. HPLC profile of the pure cs-oxidant, eluted at the retention time of 8.808 min.

10076033.052802



Fig 5. Thin layer chromatography of the pure cs-oxidant (  $R_f=0.26$  )

208250.88097001

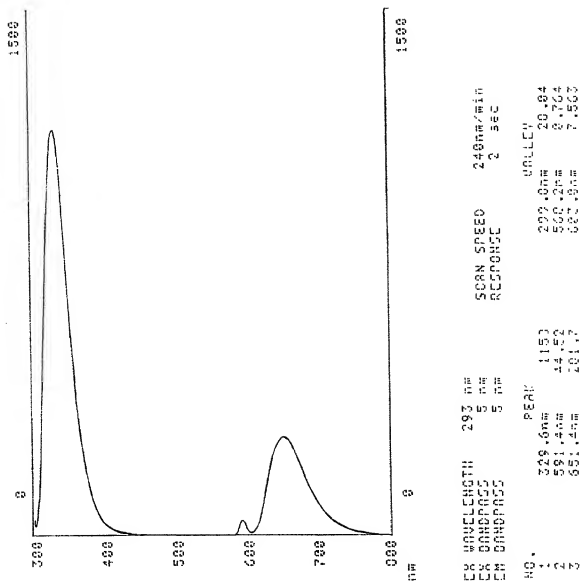


Fig 6a. Fluorescence spectroscopic profile of the cs-oxidant in methanol. The excitation was at 293 nm and the emission scanning was measured from 300 nm to 800 nm. The emission maxima were at 329.6 nm and at 651.4 nm.

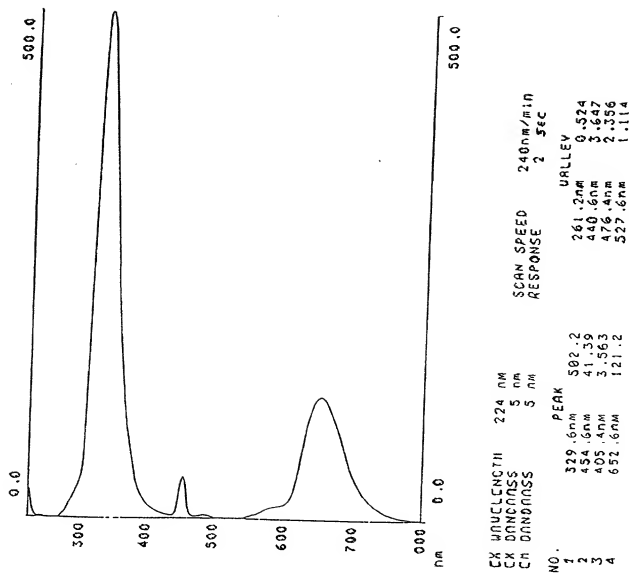


Fig 6b. Fluorescence spectroscopic profile of the CS-oxidant in methanol. The excitation was at 224 nm and the emission scanning was measured from 225 nm to 800 nm. The emission maxima were at 329.6 nm and at 652.6 nm.

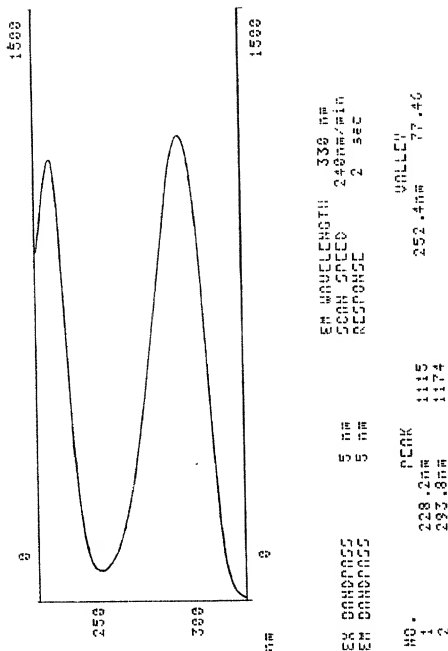


Fig 7a. Fluorescence spectroscopic profile of the cs-oxidant in methanol. The emission was at 330 nm and the excitation scanning was measured from 220 nm to 325 nm. The excitation maxima were at 228.2 nm and at 293.8 nm.



208250-EE097001

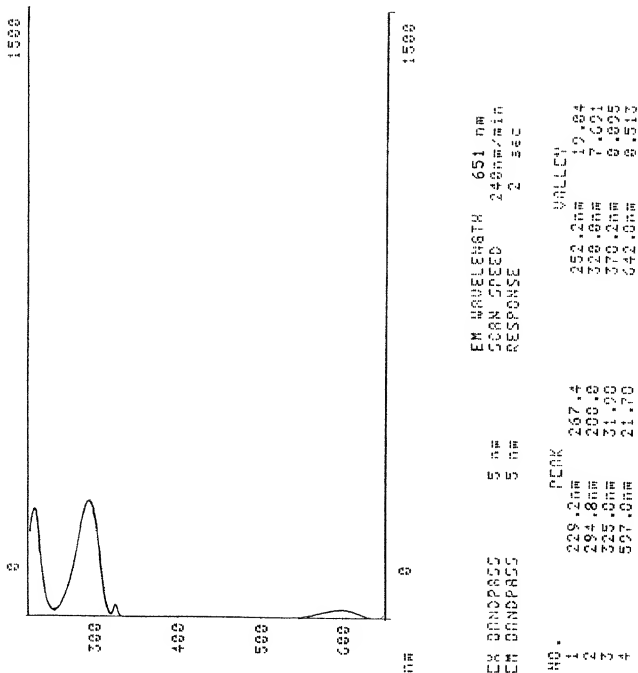


Fig 7b. Fluorescence spectroscopic profile of the cs-oxidant in methanol. The emission was at 651 nm and the excitation scanning was measured from 220 nm to 650 nm. The excitation maxima were at 229.2 nm and at 294.8 nm.



Fig 8. Crystal structure of the pure cs-oxidant

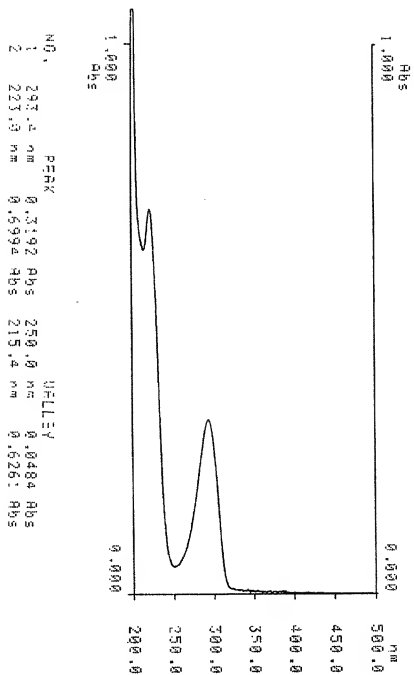


Fig 9. UV spectrophotometric profile of the cs-oxidant in methanol. It has two absorption maxima, one at 293.4 nm and another at 223.0 nm.

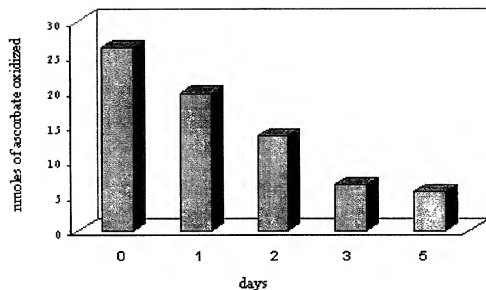


Fig.10 Stability of the solid oxidant kept at 25°C under darkness. The stability was determined by its capacity to oxidize ascorbic acid. Ascorbic acid was measured by HPLC analysis at 254 nm.

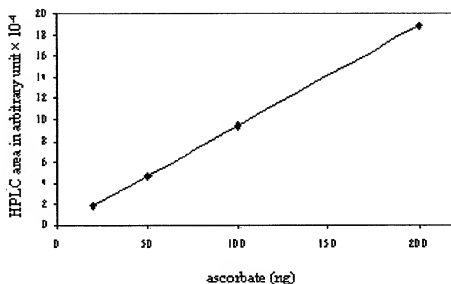


Fig.11 Standard curve of ascorbic acid based on HPLC analysis at 254 nm.

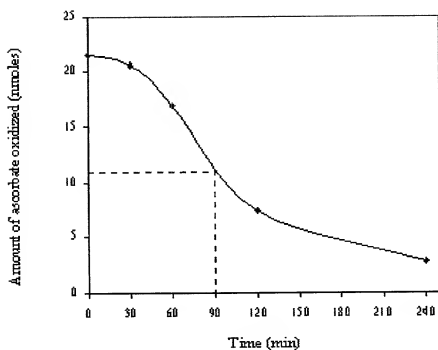


Fig.12 Stability of the cs-oxidant in 50 mM potassium phosphate buffer at 25°C measured by its potency to oxidize ascorbate as evidenced by HPLC area.

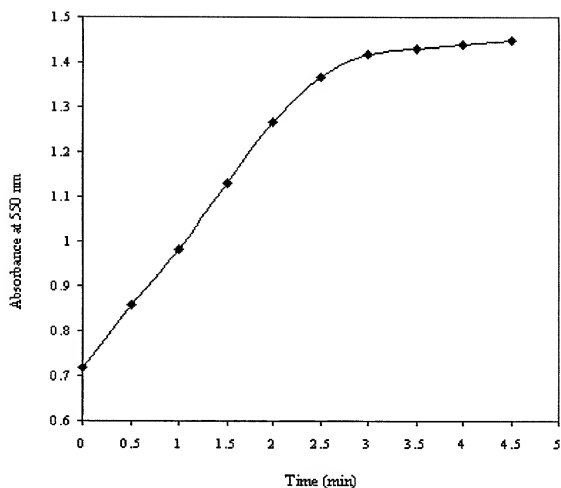


Fig.13 Quantitative reduction of ferricytochrome c by the oxidant as measured by the formation of ferrocytochrome c with time at 550 nm. The reaction was carried out in 50 mM potassium phosphate buffer, pH 7.4, keeping the final concentration of ferricytochrome c at  $100\mu\text{M}$ . One nmole of the oxidant reduced 0.71 nmoles of ferricytochrome c.

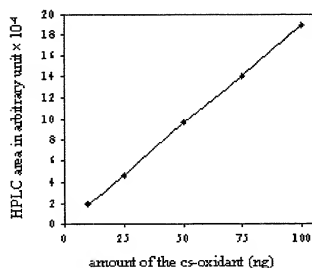


Fig 14. Standard curve of the oxidant on the basis of HPLC area at 294 nm. Different amounts of the cs-oxidant were used ranging from 10 ng to 100 ng in 20  $\mu$ l of mobile solvent.

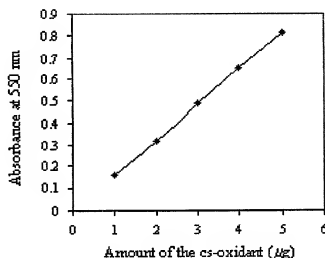


Fig 15. Standard curve of the oxidant on the basis of reduction of cytochrome c by using different amounts of the oxidant ranging from 1  $\mu$ g to 5  $\mu$ g.

208250.52802

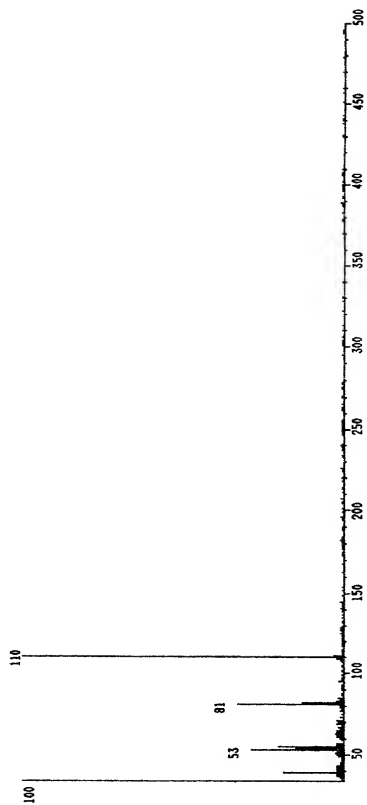


Fig 16. Mass spectrum of the pure os-oxidant



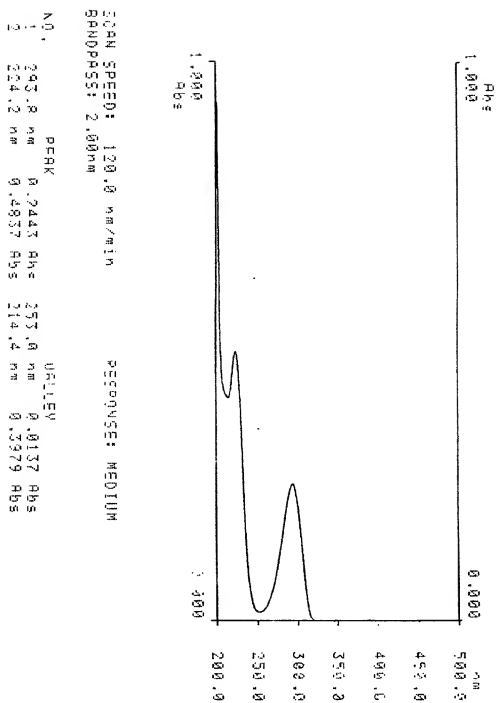


Fig 17. UV spectrophotometric profile of the hydroquinone in methanol. It has two absorption maxima, one at 293.8 nm and another at 224.2 nm.

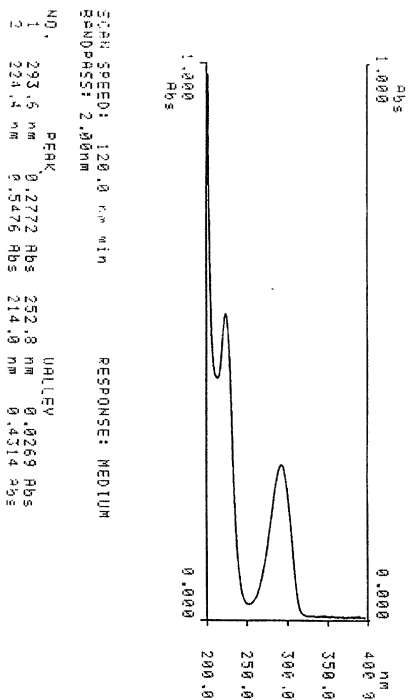


Fig 18. UV spectrophotometric profile of the cs-oxidant stored at room temperature in dark for 8 days. The two absorption maxima are at 293.6 nm and at 224.4 nm.

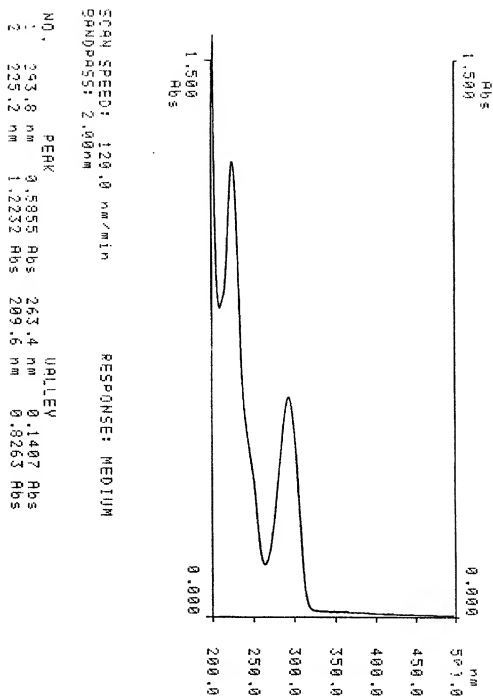


Fig 19. UV spectrophotometric profile of equimolar mixture of p-benzoquinone and hydroquinone in methanol. There is a shoulder near 242 nm ( the  $\lambda_{max}$  of p-benzoquinone).

208250.2509/001

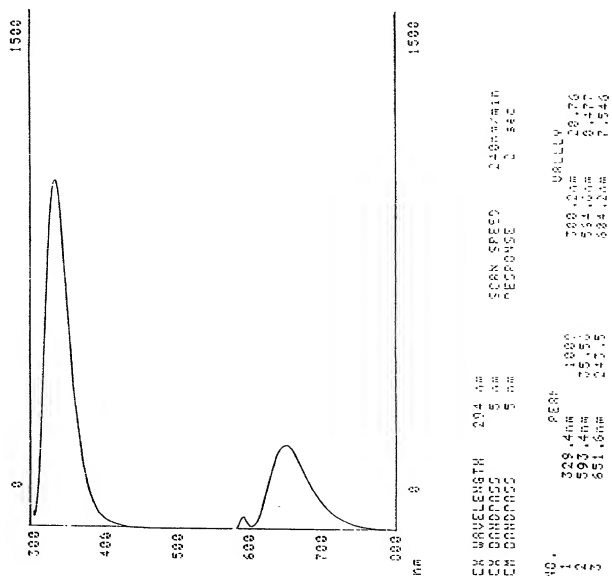


Fig 20. Fluorescence spectroscopic profile of the hydroquinone in methanol. The excitation was at 294 nm and the emission scanning was measured from 300 nm to 800 nm. The emission maxima were at 329.4 nm and at 651.6 nm.

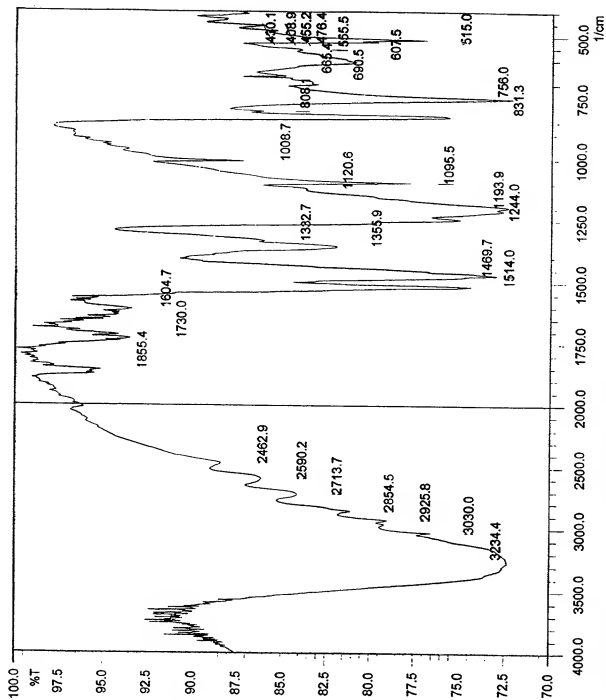


Fig 21. FTIR spectroscopic profile of the cs-oxidant.

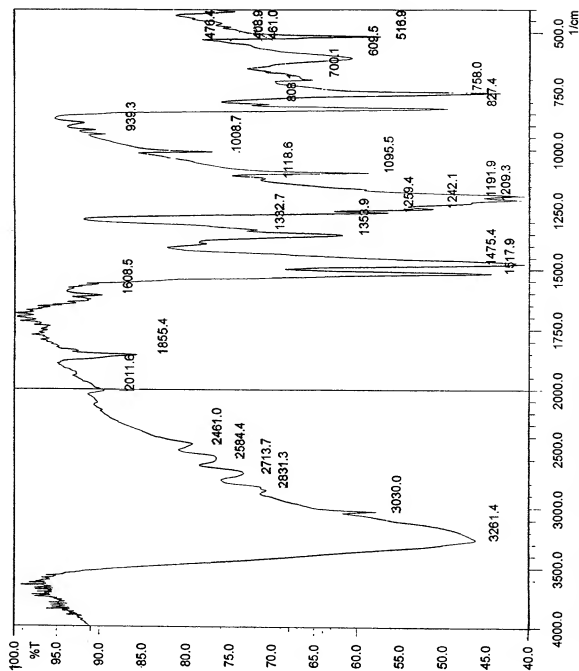


Fig 22. FTIR spectroscopic profile of hydroquinone.

10076033-052802

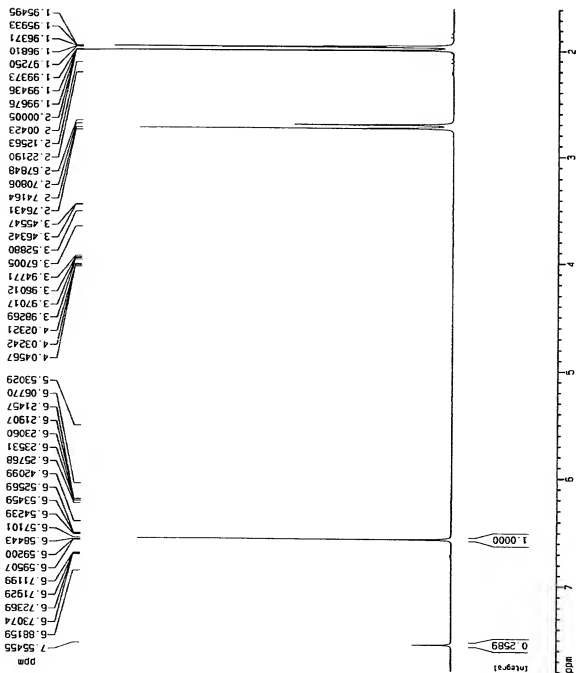


Fig 23.  $^1\text{H}$ -NMR spectroscopic profile of the ox-oxidant in  $\text{CD}_3\text{COCD}_3$

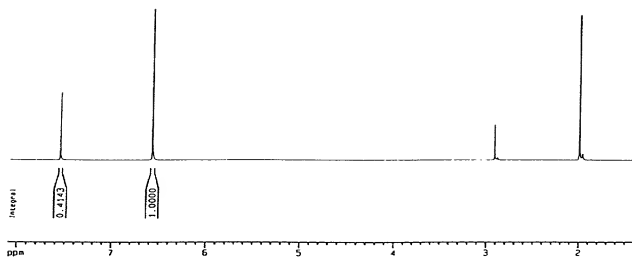


Fig 24. H-NMR spectroscopic profile of hydroquinone in  $\text{CD}_3\text{COCD}_3$

20250220/001  
10076033-052802



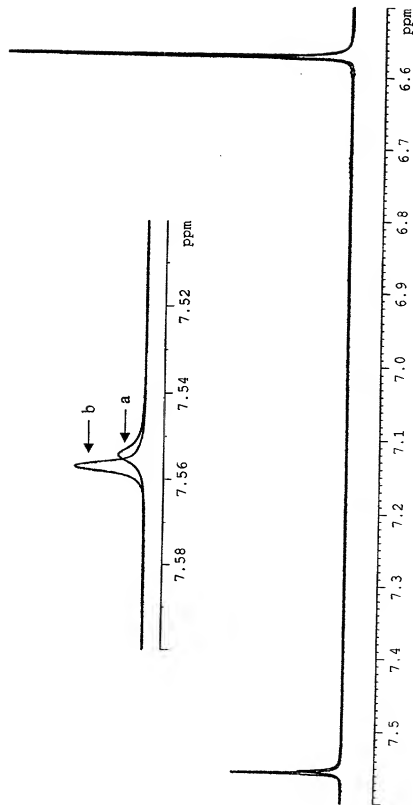


Fig 25. Comparative  $^1\text{H}$ -NMR spectroscopic profiles of (a) cs-oxidant and (b) hydroquinone.

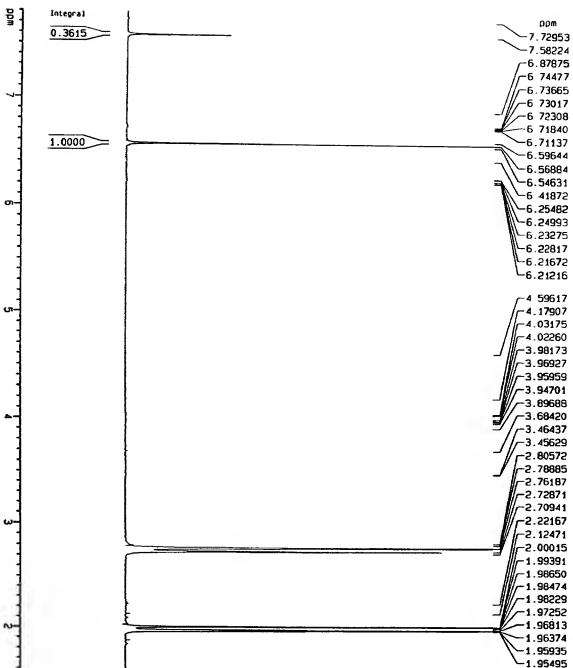


Fig 26.  $^1\text{H}$ -NMR spectroscopic profile of the ox-oxidant after reduction with sodium dithionite.

10076033.052802

10076037.052800

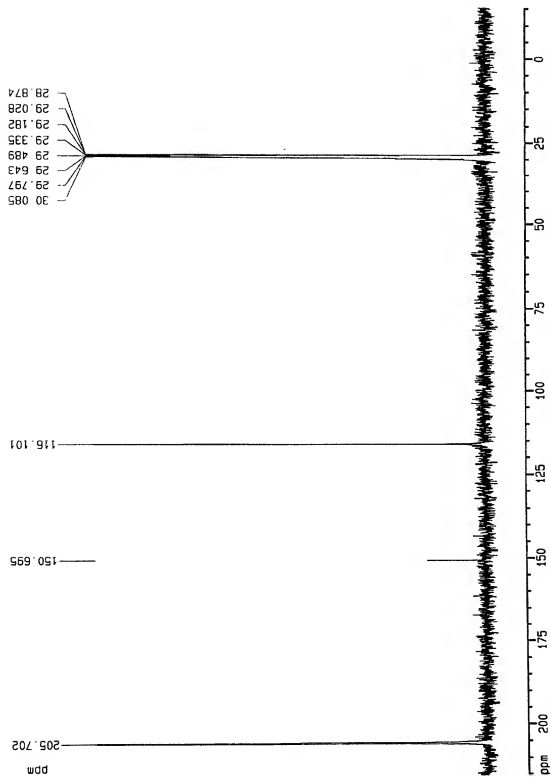


Fig 27. C-NMR spectroscopic profile of the cs-oxidant in  $\text{CD}_3\text{COCD}_3$

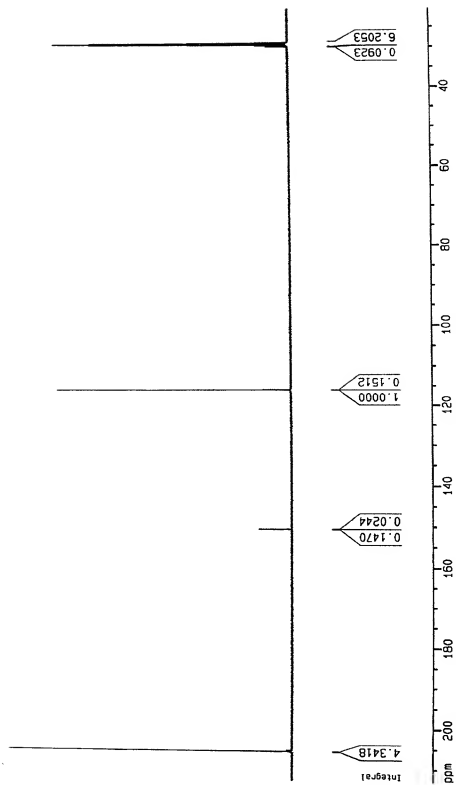


Fig 28.  $^{13}\text{C}$ -NMR spectroscopic profile of hydroquinone in  $\text{CD}_3\text{COCD}_3$

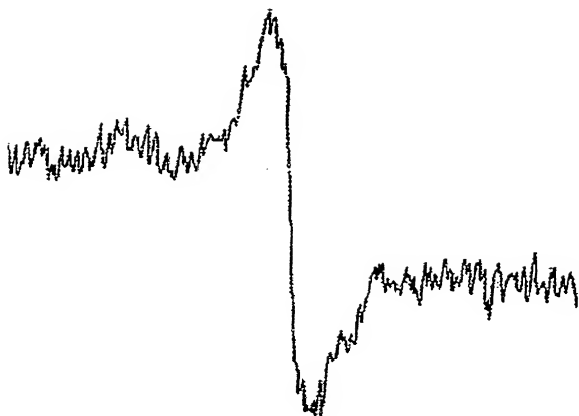


Fig 29. Room temperature ESR spectrum of CS-oxidant, freshly prepared from 100 cigarettes. The spectrum was recorded on a JES-REIX ESR spectrometer ( Tokyo, Japan). The spectral parameters were as follows: microwave frequency, 9.435 GHz; power, 2mW; field modulation width, 0.4mT; modulated frequency, 100kHz; time constant, 0.3 sec; scan rate, 2.5 mT/ sec

10076033.052802


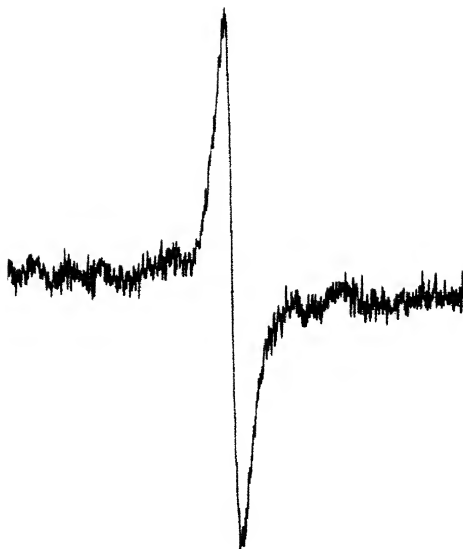


Fig 30. Room temperature ESR spectrum of aged (10 days) cs-oxidant, prepared from 400 cigarettes

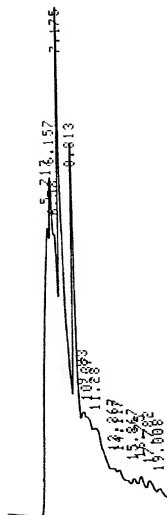

CHROMATOPAC C-6A				FILE	0	41
SAMPLE NO 0				METHOD		
REPORT NO 48						
PKNO	TIME	AREA	NK	IDNO	CONC	NAME
1	5.217	475376			19.0777	
2	6.157	317530	V		12.7431	
3	6.58	209664	V		8.4142	
4	7.175	708579	V		28.4366	
5	8.813	340583	V		13.6682*	
6	9.83	99028	V		3.9742	
7	10.37	103590	V		4.1573	
8	11.28	178509	V		7.1639	
9	13.367	24236	V		0.9727	
10	14.117	15200	V		0.61	
11	16.75	9187			0.3687	
12	17.782	10303			0.4135	
TOTAL					100	

Fig 31. HPLC profile of the whole cs solution analyzed in the silica column (Lichrospher® Si 60, Merck).

\* Indicates the retention time, area and the concentration ( 13.6682%) of the cs-oxidant

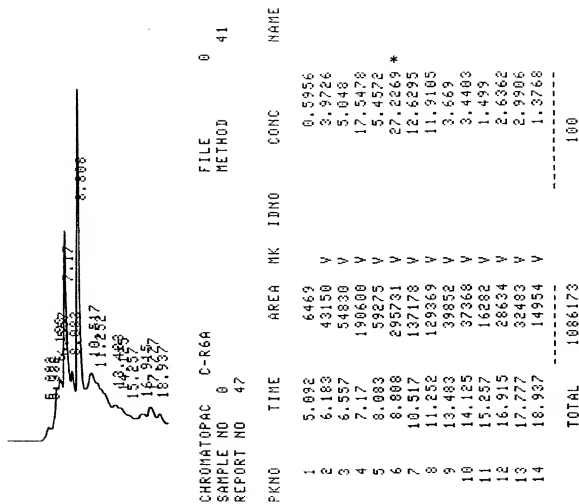


Fig 32. HPLC profile of the aqueous extract of cs solution analyzed in the silica column ( Lichrospher® Si 60, Merck).

\* Indicates the retention time, area and the concentration ( 27.2269%) of the cs-oxidant



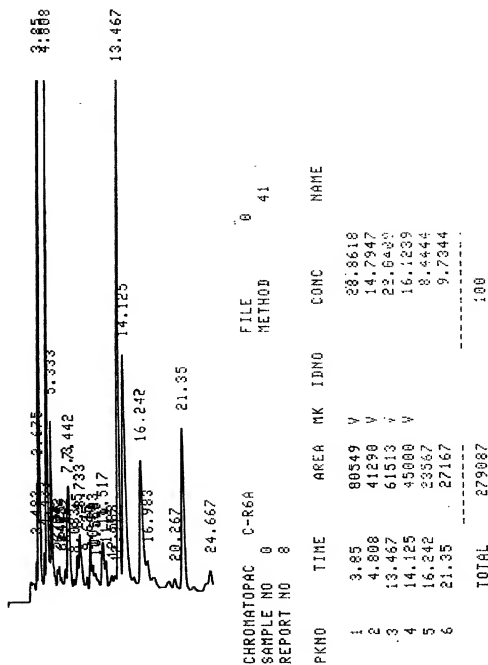


Fig 33. HPLC profile of the whole cs solution analyzed in the ODS column (Shim-pack CLC-ODS, Shimadzu). The cs-oxidant eluted at 13.467 min.

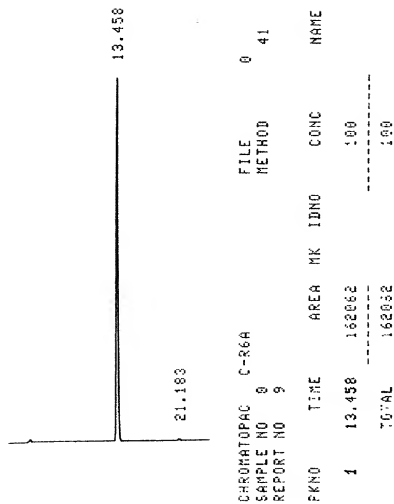


Fig 34. HPLC profile of the pure cs-oxidant, analyzed in the ODS column ( Shim-pack CLC-ODS, Shimadzu) eluted at the retention time of 13.458 min.



Fig 35a. SDS-PAGE of the guinea pig lung microsomal proteins treated with whole cs solution and the cs-oxidant. Lane 1, untreated microsomes; lane 2, microsomes treated with 50  $\mu$ l cs solution; lane 3, microsomes treated with 100  $\mu$ l cs solution; lane 4, microsomes treated with 10 $\mu$ g cs-oxidant; lane 5, microsomes treated with 20 $\mu$ g cs-oxidant.

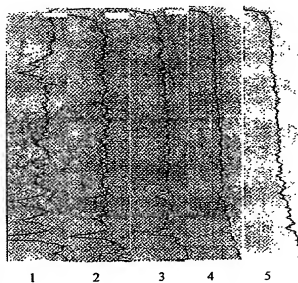


Fig 35b. Densitometric scanning of the protein bands of different lanes as in Fig 35a.